
Castellano	2	Deutsch	6
English	3	Italiano	7
Português	4	Česky	8
Français	5	Polski	9

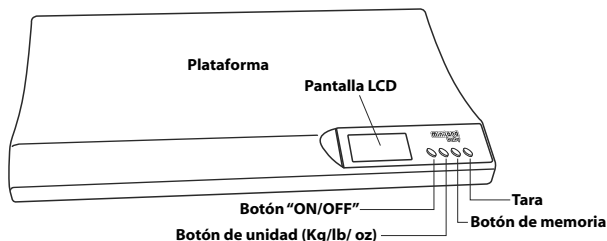Características:

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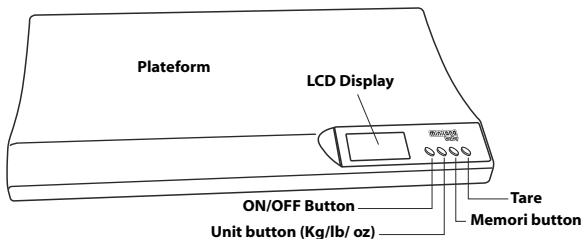
Preparations:

1. Battery installation: open the battery cover at the bottom of the scale and insert 4 brand-new AAA batteries properly.
2. Operating Environment: lay the scale onto a firm and horizontal ground.

Operation:

DO NOT PUT ANYTHING ON THE SCALE BEFORE TURNING ON. Press "ON/Off" button to turn on the scale.

Previous record appears when the scale is turned on and Do Not put the baby on the scale at this moment.

Put the baby on the scale when "0.00" appears and "(-0-)" icon blinks.

Press "Unit" button to change the measure unit (kg/lb/oz).

"HOLD" icon on the screen will blink for 5 seconds and weight is determined when it stops blinking.

To save the record, press "M" button until the icon on the screen stops blinking. New record is stored for next weight comparison.


"FFFF" appears if capacity of the scale is exceeded.

To weigh a second baby in the meantime, press "ON/Off" button and hold for 3 seconds to turn off the scale and repeat the above steps to restart.

Warning: Always be attentive to the baby during the weighing process.

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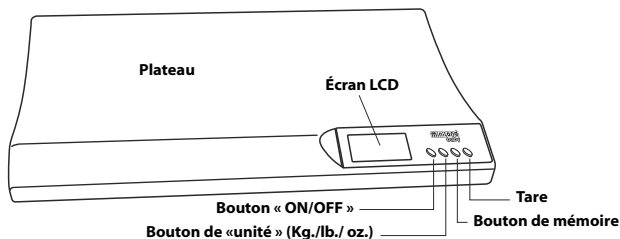We declare under our own responsibility that the product:

scaly · Miniland Baby · 89025 · Miniland, S.A.

to which this declaration refers conforms with the relevant standards or other standardising documents:

- Electromagnetic Compatibility (EMC) Directive 2004/108/EC
- EN 61000-6-1:2007 Electromagnetic compatibility (EMC). Generic standards. Immunity for residential, commercial and light-industrial environments
- EN 61000-6-3:2007 Electromagnetic compatibility (EMC). Generic standards. Emission standard for residential, commercial and light-industrial environments
- Waste Electrical and Electronic Equipment (WEEE) Directive 2002/96/EC

Miniland, S.A.
Chief Executive Officer
José Verdú Francés
Onil, Alicante, Spain.
31st March 2011

Miniland S.A. P.Ind. La Marjal I C/ La Patronal s/nº
03430 ONIL (Alicante) SPAIN

☎ Technical support: +34 966 557 775

Fax +34 965 565 454

www.minilandbaby.com · miniland@miniland.es

© **Miniland, S.A. 2009**

Fabricado en China. Conservar los datos.
Made in China. Keep this information.
Fabricado em China. Conservar dados.
Fabbricato in China. Conservare i dati.
Fabriqué en China. Conserver ces informations.
Hergestellt in China. Bewahren Sie diese Angaben auf
Wyprodukowano w Chinach. Zachowaj instrukcję.
Vyrobeno v Číně. Tyto informace si prosím uchovejte.
Сделано в Китае. Сохраните эту информацию.

Recogida selectiva. Cuida el medio ambiente.
Selective pick-up. Protect the environment.
Tri sélectif. Prend soin de l'environnement
Bitte gesondert entsorgen. Schützen Sie die Umwelt.
Raccolta differenziata. Rispetta l'ambiente.
Recolha selectiva. Conserve o Meio Ambiente.
Chroń środowisko! Segreguj śmieci!
Tříděný odpad. Chráníme životní prostředí.
Раздельный сбор мусора. Берегите окружающую среду.

89025 scaly